


Application reference

The Solaris® gold bath can be used together with the Solaris® activator in a Wieland AGC Micro Plus unit.


Care must be taken to ensure the mixing ratio of 100 ml of gold bath to 5 ml of activator is maintained.

The amount of gold bath required for precipitation can be determined by consulting Wieland's table for the Wieland AGC Micro Plus.

All other device parameter settings can be determined from the corresponding original instruction manual of the unit.

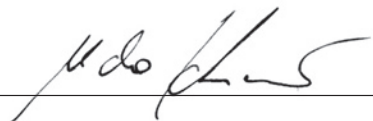

Dr. Udo Schusser
Director, Research and Development

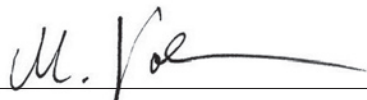

Dr. Markus Vollmann
Director, Application Technology