

Executive[®]

Aesthetic and harmonic.
Setting the tone in prosthodontics.



DENTSPLY

Setting the tone in prosthodontics

To realize this ambitious goal, we have mobilized all our resources: intensive research, valuable raw materials, up-to-date production methods, and rigour quality control.

Executive®

Aesthetic and harmonic.








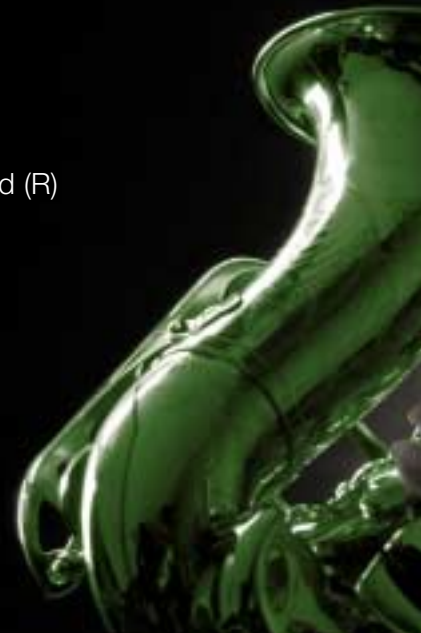
Aesthetics and rationalism in harmony

The Executive® denture teeth made of INPEN® are a perfect answer to demands for aesthetic sophistication and economy alike. Their special shade design gives these teeth a natural opalescence and vivid appearance. The posteriors low cusp inclination (20° – 22°) allow easy set-up. A good central occlusion is ensured right from the start. Masticatory balance is easy to achieve, minimizing the need for adjustments at the laboratory and at chairside – and forming the basis for rational and efficient procedures.


Transparent mould coding

Executive® anteriors are subdivided into five mould classes based on the type-specific theory of harmony, making the right choice of teeth much easier.

1.  Square (S)
2.  Square Ovoid (SO)
3.  Square Tapering (ST)
4.  Tapering (T)
5.  Round/Ovoid (R)



INPEN®-Structure:
Homogenous three-dimensional structure and connection of primary and secondary molecular chains

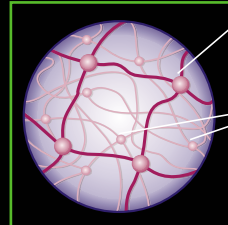

Connection of primary chains

Connection of secondary chains



Economical range of moulds

Twenty-one moulds are available for the maxilla and eleven for the mandible. This range is supplemented by five posterior moulds. All Executive® teeth are available in V-shades (A1 – D4).

INPEN® – Quality Acrylic

The high-quality INPEN® material (interpenetrated polymer network) is the chemical basis of the DENTSPLY denture teeth. Its highly connected copolymers give this material its superior properties:

- Excellent abrasion resistance
- Minimal plaque affinity
- Excellent colour stability
- Excellent craquelé resistance
- Maximum bond strength
- Added reliability for the laboratory, dentist and patient



The result:
high-quality and long-lasting
denture teeth.



Executive®

Upper anteriors

Available in V-shades (A1–D4)






Mould	Image	Width of the complete set in mm	Length of the centrals without neck in mm	Width of the centrals in mm
	Square			
S1		44.0	9.3	8.1
S2		45.5	10.6	8.5
S3		49.0	10.6	9.2
S4		48.0	11.1	8.8
S5		50.5	11.6	9.4





Mould	Image	Width of the complete set in mm	Length of the centrals without neck in mm	Width of the centrals in mm
	Square Ovoid			
SO0		40.5	9.2	7.5
SO1		45.0	10.0	8.2
SO2		46.0	11.5	8.6
	Square Tapering			
ST1		44.0	10.5	8.1
ST2		45.0	10.8	8.5
ST3		46.5	10.9	8.8
ST6		54.5	12.4	10.2

Executive®

Upper anteriors

Available in V-shades (A1–D4)

Mould	Tapering	Width of the complete set in mm	Length of the centrals without neck in mm	Width of the centrals in mm
T1		45.5	9.6	8.8
T2		43.5	10.3	8.0
T3		48.0	10.4	8.7
T4		46.5	10.5	8.5
T7		53.5	13.3	10.0

Mould	Round/Ovoid	Width of the complete set in mm	Length of the centrals without neck in mm	Width of the centrals in mm
R1		43.5	10.3	8.3
R2		45.5	11.0	8.4
R25		46.5	10.9	8.3
R3		48.5	11.9	9.2

Executive®

Lower anteriors

Available in V-shades (A1–D4)

Mould

Width of the complete set in mm
Length of the centrals without neck in mm
Width of the centrals in mm

Mould	Width of the complete set in mm	Length of the centrals without neck in mm	Width of the centrals in mm
L1	35.5	8.8	5.4
L2	32.0	8.9	5.0
L3	32.0	9.1	4.9
L4	33.5	9.1	5.2
L5	36.0	9.3	5.2
L6	33.5	9.4	5.0
L7	38.5	9.7	5.7
L8	35.5	9.8	5.4

Mould

Width of the complete set in mm
Length of the centrals without neck in mm
Width of the centrals in mm

Mould	Width of the complete set in mm	Length of the centrals without neck in mm	Width of the centrals in mm
L9	35.5	10.3	5.5
L12	38.5	11.9	5.9
L13	37.0	12.8	5.8

Executive®

Posteriors

Available in V-shades (A1–D4)

Mould

U31



L31

U32



L32

U33



L33

U34




L34

U36



L36

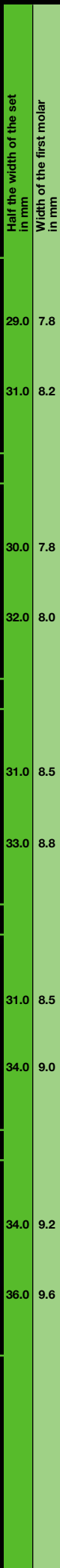

Exact shade matching



With the original Executive® shade guide

in V-shades (A1–D4)

Order No. D 01 027 215

Combination table

Executive® Upper anteriors	Executive® Lower anteriors	Executive® Posteriors
Square		
S1	L3	33
S2	L4	33
S3	L5	34/36
S4	L5	33
NEW: S5	L7	34/36
Square Ovoid		
NEW: SO0	L2	31
SO1	L3	33
SO2	L7, L8, L9	34/36
Square Tapering		
ST1	L2	31/32
ST2	L3	31/32
ST3	L1	33
NEW: ST6	NEW: L12, L13	34/36
Tapering		
T1	L4	33
T2	L2	31/32
T3	L1, L8	33/36
T4	L3	33
NEW: T7	NEW: L12, L13	36
Round/Ovoid		
R1	L6	32
R2	L4	32
NEW: R25	L1, L8	32
R3	L5	34/36

